



OB/GYN

Electrode	Surgery	Function	Key concept	Standard accessories	Option accessories
E01	Cervical erosioin / CIN	Coagulation	The more wider diameter , The deeper coagulated		
E02					
E03					
E04					
E06					
E07	Cervical erosion /CIN	Coagulation & Hemostatis	Hemostatic cutting and surface coagulation		
E09					
F04	Cervical erosion /CIN	Cutting in cone shape (Conization)			
J03	Cervial erosion /CIN	pure cutting			
K02	Cervical erosion /CIN	Tumor removal			
K04		Biopsy and removal			
L02					
N002					
O01					
G03	Round LETTZ (LEEP)		Loop Electrosurgical Excision Procedure		


* **F 04 (Cornization electrode for cervical erosion , CIN)**

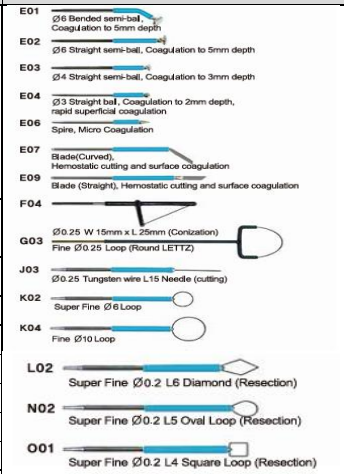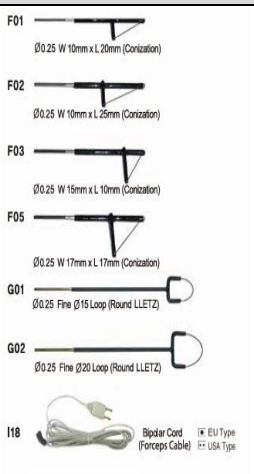

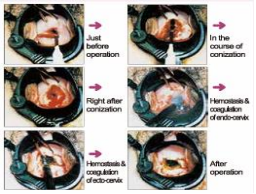


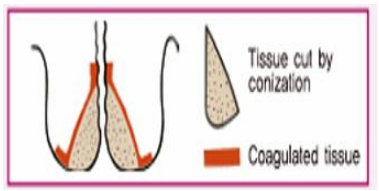

ENT

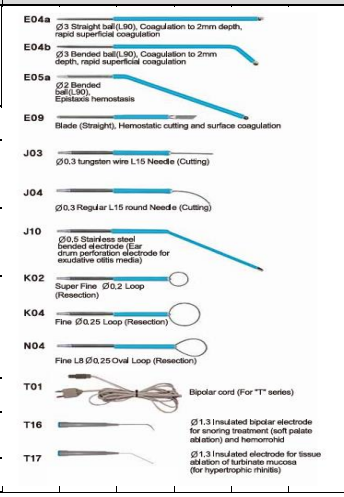
Electrode	Surgery	Function	Key concept	Standard accessories	Option accessories
E04a		Coagulation	The more wider diameter, the deeper coagulated		
E04b					
E05a					
E09		Coagulation & Hemostatis	Hemostatic cutting and surface coagulation		
J03		Pure cutting			
J04					
J10	Ear drum perforation Exudative ottis media				
K02					
K04					
N04					
T16	Snoring				
T17	Inferior Turbinate Hypertrophy				
T01			For T16, T17 connection (Bipolar)		

T 16 (Snoring)

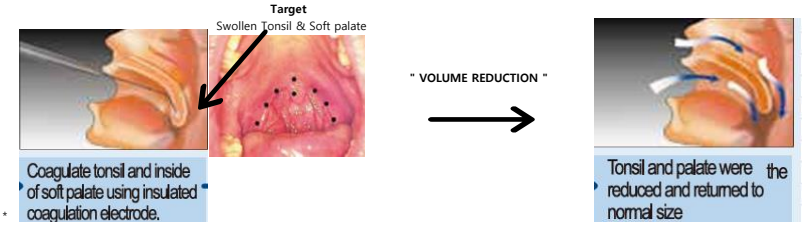

*

T17 Inferior Turbinate Hypertrophy

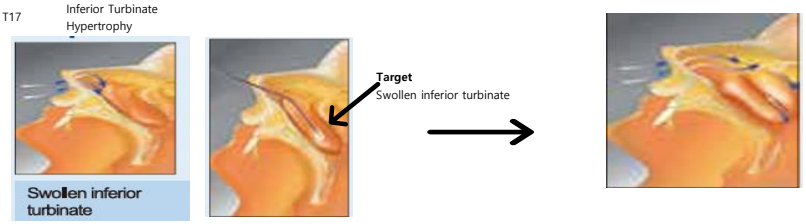

Dermatology & Urology

Electrode	Surgery	Function	Key concept	Standard accessories	Option accessories
E04		Coagulation	The more wider diameter, the deeper coagulated	E04 Ø3 Straight ball(Hemostatic coagulation) E05 Ø2 Straight ball(Hemostatic coagulation) E06 Spire(Micro coagulation) E09 Blade (Straight), Hemostatic cutting and surface coagulation E10b Hemostatic cutting electrode for blepharoplasty (15mm) J03 Regular L15 Needle(Incision, Vaporization) J06 Super fine L5 bended Needle J11 Micro cutting (2mm) J12 Micro cutting (3mm) K01 Super fine Ø4 Loop K02 Super fine Ø6 Loop	BS2510 Ø0.15 wire(Hair root removal) BS3512 Ø0.15 wire(Hair root removal) BS5025 Ø0.15 wire(Hair root removal) BL4015 Ø0.18 wire(Hair root removal) BL4515 Ø0.18 wire(Hair root removal) BL5025 Ø0.18 wire(Hair root removal) E10a Hemostatic cutting electrode for blepharoplasty (10mm) E10c Hemostatic cutting electrode for blepharoplasty (20mm) I18 Bipolar cord (Forceps Cable) EU Type USA Type
E05a					
E06					
E09		Coagulation & Hemostatis			
E10a					
E10C					
E10b	Circumsizing	Pure cutting			
J03					
J06					
J11		Micro cutting			
J12					
K01	Polyp, wart , freckle removal	Pure cutting			
K02					
BS2510	Hair root removal	Micro coagulation			
BS35112					
BS5025					
BL4015					
BL4515					
BL5025					
I18			Bipolar forcep connection		

"Polyp Removal by resecting"




" Freckle removal by vaporizing"



" Circumsizing by purcutting "



"Hair root removal"

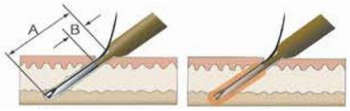

GS (Proctology)

Electrode	Surgery	Function	Key concept	Standard accessories	Option accessories
E03b	Hemorrhoid	Coagulation	The more wider diameter, the deeper coagulated	E03b Ø4 Bended semi-ball(L90), Coagulation to 3mm depth E04a Ø3 Straight ball(L90), Coagulation to 2mm depth, rapid superficial coagulation E04b Ø3 Bended ball(L90), Coagulation to 2mm depth, rapid superficial coagulation E06 Spire, Micro coagulation E07 Blade(Curved), Hemostatic cutting and surface coagulation E09 Blade (Straight), Hemostatic cutting and surface coagulation J03 Ø0.3 Regular L15 Needle J04 Ø0.3 Regular L15 Round Needle	I18 Bipolar Cord (Forceps Cable) EU Type USA Type T01 Bipolar cord (For "T" series) T16 Ø1.3 Insulated bipolar electrode for snoring (soft palate cogulation) and hemorrohid
E04a					
E04b					
E06		Micro coagulation			
E07	Hemorrhoid	Coagulation & Hemostatis	Hemostatic cutting and surface coagulation		
E09					
J03		Pure cutting			
J04					
T16	Hemorrhoid		Volume reduction of Hemorrhoid		
I18			Bipolar forcep connection		
T01			For T16 connection (Bipolar)		

Hemorrhoid surgery

